小米體重計2 使用說明書	-1
Mi Smart Scale 2 User Manual	-6
Mi Smart Scale 2 Manual de usuario	-11
Mi Smart Scale 2 Mode d'emploi	-16
Mi Smart Scale 2 Manuale utente	-21
Mi Intelligente Waage 2 Benutzerhandbuch	-26
Весы Mi Smart Scale 2 Руководство пользователя	-31
Mi 스마트 체중계 2 사용 설명서	-36
Inteligenta waga Mi 2 Instrukcja obsługi	-41
Mi スマート体重計 2 取扱説明書	-46

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
Read this manual carefully before use and retain it for future reference.

Installing the Batteries

Remove the sticker and open the battery cover, then insert three 1.5V AAA batteries as indicated by the illustration.

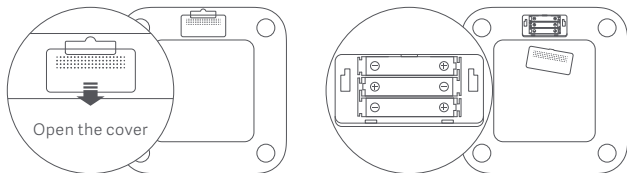Installing the app

Scan the QR code to download and install the app. Add the Mi Smart Scale 2 to the app first before you start using it.

Pairing the Scale

Open the app and sign in using your Mi account, then pair the scale with your account and enter your personal data according to the on-screen instructions. After pairing the scale, you can check the weighing records and weight analysis via the app.

Personal data is used to analyze your weight. To ensure the accuracy of the results, please enter accurate personal data.

Weighing People

Place the scale on a hard, flat floor.

Follow the orientation and position shown in the figure, step on the scale and stand steadily on both feet in the middle of the scale.

Once the indicated weight blinks, the weighing has been completed and the weight analysis information can be viewed on your phone.

The weight analysis information provided by this scale is for reference only and does not constitute medical advice.

Weighing Items

Place the scale on a hard, flat surface.

Press on the scale to turn it on. Wait until "0.00" is displayed and place the item you want to weigh onto the scale.


Moving or bumping the scale may affect the accuracy, in which case you should weigh the item again to make sure the measurement is accurate.

Alert Indications

When weighing, a blinking "L" indicates a low battery level, and you should replace the batteries as soon as possible.

When weighing, a blinking "OL" indicates the scale is overloaded beyond its maximum capacity of 150 kg.

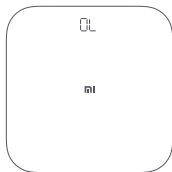Selecting a Weighing Unit

The scale calculates weight in kilograms (kg) by default. You can change the scale's weighing unit via the device settings in the app.

The scale can display weight in kilograms (kg), and pounds (lbs).

Precautions

1. To prevent slipping, do not place the scale on a slippery floor, and do not use it when your feet or the scale are wet.
2. Always stand in the middle of the scale, to prevent it from accidentally flipping over while weighing yourself.
3. Only use the scale in a dry environment, and do not immerse it in water.
4. Wipe the surface of the scale with a damp cloth and mild detergent to clean it. Do not use any organic solvents or other chemicals.
5. After the batteries are depleted, immediately remove them to avoid potential leaking and corrosion.
6. Do not replace a battery with an incorrect type. Disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery, that can result in an explosion. Leaving a battery in an extremely high temperature surrounding environment that can result in an explosion or the leakage of flammable liquid or gas. A battery subjected to extremely low air pressure that may result in an explosion or the leakage of flammable liquid or gas.
7. This scale is not recommended for persons with a pacemaker or other implanted medical devices. The accuracy of the body composition data may vary when used by professional athletes, pregnant women, and certain other individuals with special conditions.
8. You can go to app to experience the function of One-foot-stand with Eyes Closed Measurement. When you are feeling off balance, please be sure to open your eyes and stand on your feet firmly to prevent from accidentally flipping over and causing a fall. Anyone including children, elderly, expectant mother, and others with physical or balance disorders should use this function with caution.

Disposal and Recycling Information

The following symbol signifies that this product has to be disposed separately from ordinary household wastes at its end of life. Please kindly be aware that this is your responsibility to dispose the electronic equipment at recycling centers so as to help conserve the natural resources.

Certifications and Safety Approvals

Hereby, Anhui Huami Information Technology Co., Ltd. declares that the radio equipment type XMTZC04HM is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://www.mi.com/global/service/support/declaration.html>

The product is in compliance with UK PSTI regulations, the full text of declaration of conformity is available at the following internet address:
<https://www.mi.com/global/service/support/declaration.html>

Specifications

Name: Mi Smart Scale 2

Model: XMTZC04HM

Dimensions: 280 × 280 × 22 mm

Net Weight: 1.2 kg

Weighing Range: 100 g - 150 kg

Weighing Units: kg / lbs

Increments: 0.1kg/lbs

Compatible with: Android 4.4 or iOS 10.0 (iPhone 5S and later) and above, and Bluetooth 4.0 and above

Casing Material: Tempered glass, engineering plastic

Wireless Connectivity: Bluetooth 5.0

Operating Temperature: 0°C to 40°C

Power Source: 1.5 V AAA batteries ×3

Frequency: 2402-2480MHz

Max. Output: 4 dBm

Mi Smart Scale 2

Model: XMTZC04HM

Manufactured for: Xiaomi Communications Co. Ltd.

Manufactured by: Anhui Huami Information Technology Co., Ltd.

Address: 7/F, Building B2, Huami Global Innovation Center, No. 900, Wangjiang West Road, High-tech Zone, Hefei City, China (Anhui) Pilot Free Trade Zone

For further information, please go to www.mi.com

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